

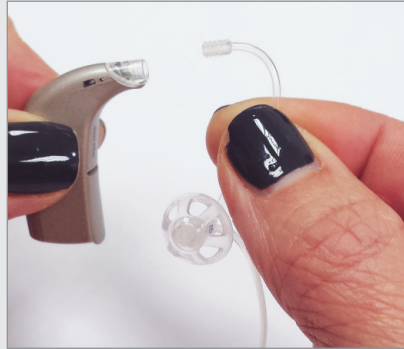


Cleaning the Corda miniFit (thin tube)




The following pictures illustrate how to clean a hearing aid with a thin tube.

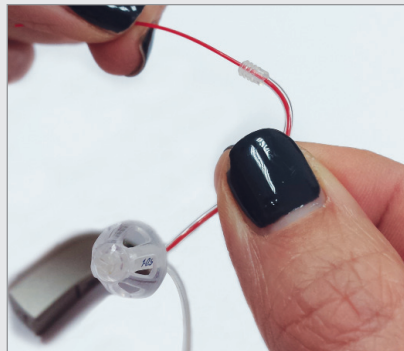


1. Detach the corda tube by pulling the tube away from the hearing instrument.

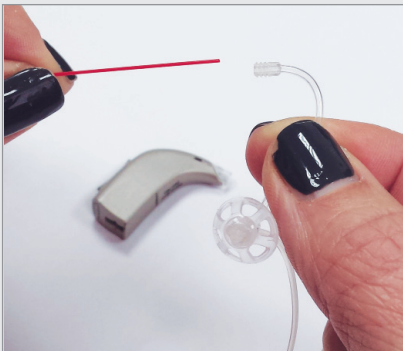

2. Insert the cleaning tool into the top part of the tube and feed it all the way through to the dome.

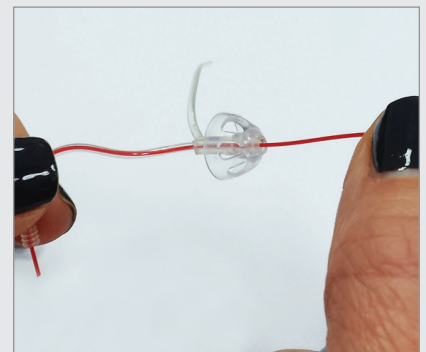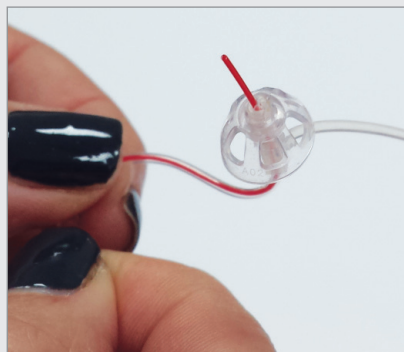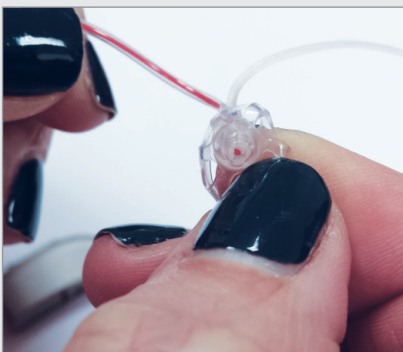


3. Carefully push the cleaning tool through the dome past the vent, without removing the dome.

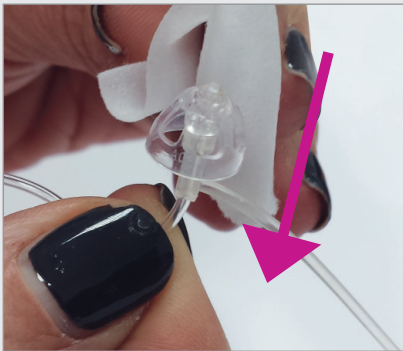

*Always clean the dome
from tip to base*

4. To clean the dome, wipe the dome gently with a wet-wipe or damp cloth to remove any excess cerumen (wax).

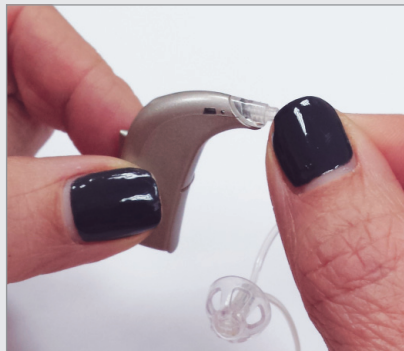

5. Re-attach the corda tube by pushing it back on to the hearing instrument with a simple click fit action.

PLEASE NOTE:

1. If the dome is removed, always replace with a NEW dome.
Never reattach a used dome back onto the tube.